medivon.

User manual
MESSAGE MAT

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1. Setup and operation

1. Secure the cushion in place by wrapping the back elastic straps over the car seat or chair.
2. Connect the adapter cable to the corresponding cable at the side of the cushion.

3. Plug the home adapter into an outlet or the car adapter into your car's power receptacle.
4. Turn on the device using the controller.
5. When finished place the controller in the pouch at the side of the seat cushion and unplug the adapter from the outlet or receptacle.

- Do not setup or use the massager in a bathroom or similiary wet/damp areas as such an environment may reduce the lifespan of the massager mechanism.
- Ensure the massager is not exposed to direct sunlight or placed near sources of heat.


Note: This product produces a warming heat. It is not intended to be used as a therapeutic heating pad. Please allow a few minutes to reach maximum temperature. The cushion reaches a maximum temperature of 51°C (124°F) after approx. 7 minutes.

2. Safety Instructions

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- Do not wet, do not use pins, never remove cover.
- This items is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff.
- Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 60 minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.

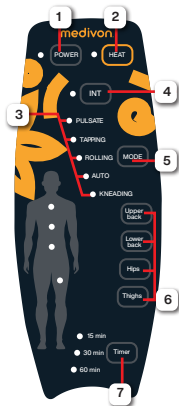3. Controller instructions

1. Power:

Press to turn the massager on or off.

2. Heat:

Press to activate heat.

3. Mode Indicators; Mode Descriptions:

Pulsate: Vibration increases and decreases at regular intervals.

Tapping: Short bursts of vibration.

Rolling: Moves through each area, from upper back to thighs, and then repeats at short intervals.

Auto: Cycles through various massage patterns.

Kneading: Constant vibration.

4. Intensity:

Press to choose vibration intensity (4 levels optional).

5. MODE:

Press to change the mode.

6. Massage Regions:

Each massage region may be turned off or on as desired, for any combination of regions and while any of the five modes are active.

7. Timer Button:

Press to set the massage / heat time at 15/30/60 minutes.

4. Care and Cleaning Instructions

- Do not wash in washing machine or immerse the product or the adapter in water.
- Always switch off before cleaning and remove the adapter from the socket with clean, dry hands.
- The massager should be wiped with a slightly moist cloth and dried thoroughly. Do not wipe the product with alcohol or harsh chemicals such as thinner or bleach.
- Do not overstretch the elastic back straps.
- Keep away from heat, fire or exposure to direct sunlight.
- Store the massager in a cool, dry place out of reach of children.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health damage or environmental pollution during further disposal or recycling of the equipment.

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