
medivon.

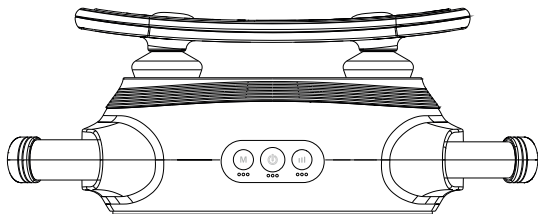FITNESS VIBRATION PLATFORM

Home Relax Therapy Series

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1. Device operation

Power button: Press and hold the ON/OFF button for 3 seconds to turn on the device, Press again to turn off the device.

Mode button: Press the mode button (M) to cycle through the type of massage to suit your individual needs. After selecting the mode, the appropriate indicator lights up.

Message intensity indicator: Press the button to rotate the message intensity according to the user's personal preferences. After selecting a level, the appropriate indicator lights up.

Landing:

1. Connect the charging cable to a DC5V 1A or DC5V 2A power adapter (usually you can use an adapter that is used to charge your phone).
2. Short press the power button to check the battery capacity when the power is off.
3. The power button and message intensity indicator will flash while charging. When the device is fully charged, the indicator light will stop flashing.

NOTE: The massager should not be used while charging.

2. Warnings

For your safety and to avoid property damage, please read this manual carefully and use as directed. Improper use may result in serious injury or death and/or property damage. Do not use this device other than as indicated in the instructions:

- Only use the charger that came with this device at the time of purchase.
- Use of an unspecified charger or adapter may result in overheating, electric shock, or fire.
- For your safety, only use the accessories and attachments provided with this massager at the time of purchase.
- DO NOT use the product while driving.
- DO NOT use the massager with wet hands.
- DO NOT insert any objects into the charger input socket. This may damage the product or cause electric shock.
- Keep the product away from heating appliances such as oven, microwave, hot cookware, pressure cooker. This may result in battery degradation and/or product failure.
- DO NOT use heavy objects or excessive pressure on the product.

- DO NOT disassemble or drill into the product and/or batteries. This may cause product failure, electric shock or fire.
- Keep fingers, hair or any other part of your body away from the shaft and back of the applicator as pinching may occur.
- Use only on a dry body surface, on clothing, pressing gently and moving over the body for approximately 60 seconds per region.
- Use the massager only on soft body tissue as desired, stop using immediately if you experience pain or discomfort. Do not use on the head or on hard or bony areas of the body.
- Bruises may occur regardless of control setting or pressure applied.
- Inspect treated areas frequently and discontinue use immediately at the first sign of pain or discomfort.
- The muscle massager is intended for external use on body parts only.

If you have any of the following conditions, consult your doctor before use:

- Pregnant women, children or elderly people
- People suffering from skin infections or having skin redness and burns who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors etc.
- People who are undergoing treatment or experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

3. Cleaning and care

Before cleaning, turn off the product and unplug the power cord.

Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with corrosive detergents such as benzene, thinners. Store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environments, otherwise it may cause colors to fade or damage the product.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.
The device was placed on the market after 13.08.2005.

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The manufacturer declares under his sole responsibility that the Medivon Force device complies with Directive 2014/30/UE of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Force