
medivon.
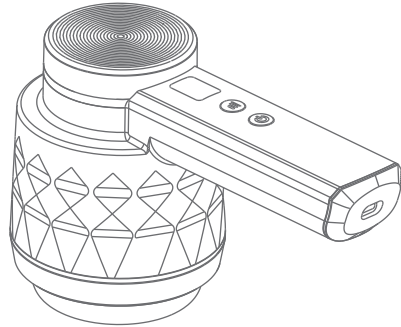
## Body Massager

Deep relaxation

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## 1. Operation and application

Turn on the massager by pressing the power button (3) for 2 seconds; the device will start and enter working mode.

In working mode, a short press of the power button (3) allows you to cycle between working modes; There are five different modes available. A short press on the adjustment button (4) allows you to adjust the massage strength; 9 intensity levels are available.

To turn off the device, hold down the power button (3) for 2 seconds. The massager is powered by a built-in lithium battery. If the battery is exhausted or low, please charge the device before next use. Use the included cable for charging. Check that the power connection is correct and that the cable is properly connected to the device (2). Voltage instability may occur during charging; Do not use the product while charging. The charging process takes approximately 4 hours. After charging, unplug the power.

Replacing massage heads:

As shown in the picture (6), turn counterclockwise to unscrew and remove the massage head mesh cover (1). After replacing the head, screw on the mesh cover clockwise. The heads can be replaced by pulling them upwards.

## 2. Dangers and warnings

People with medical implants, undergoing treatment, heart disease, high blood pressure, acute illness, pregnancy, osteoporosis, spine fractures, skin damage, fever above 38°C and people who are insensitive to heat should consult a doctor before use. Minors and people with intellectual disabilities may use the product only under supervision. The product is intended for use by one person, do not use it with others. After moving the product from a cold to a warm place, wait an hour before using it. Recommended use time is 10 minutes. Do not handle the product with wet hands and do not use it outdoors. Only use a power supply with the appropriate voltage. Do not smoke while using the product. Please use a qualified USB adapter for charging. Do not use in a humid environment. Do not use the product immediately after meals or simultaneously with other therapeutic devices. If there are any abnormalities during use, discontinue use immediately and consult a doctor. If the power cord is damaged, contact your distributor for replacement. When charging is complete or when not in use, unplug the plug from the outlet. Do not use the product while standing on it or damage the power cord. Do not use on an electric heating blanket or space heater. Do not disassemble the product unless you are an authorized service technician. In the event of a power outage, immediately turn off the product and unplug it from the power supply. Do not use the product while sleeping, intoxicated or feeling unwell. If there are any abnormalities (smoke, unpleasant odor, etc.), discontinue use immediately. Do not clean the product with gasoline, solvents or insecticides. In case of water spill, do not use the product. Do not use during a lightning storm. Before disposing of the product, please remove the battery. Keep these tips for safe and effective use of the product.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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This device was placed on the market after 13.08.2005.

The manufacturer hereby declares that the Medivon Lift device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska

Model: Medivon Lift