
medivon.
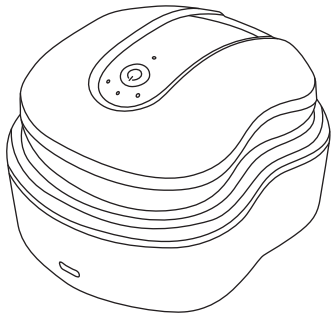
## User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK  
SL | EE | LT | LV | HU | RO | RU | NB | FI  
EL | HR | SV | DA | NL | TR | BG | AR

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## Device Manual

The device is ready for immediate use after purchase. Make sure the device is charged, then start the massager and select the appropriate operating mode according to your preferences. When the massager is turned off, short press the button to turn it on. The first slow speed mode will start. Press the button to cyclically change the operating modes or massage direction and speed.

If the device is discharged, the LED will flash red. Connect the device to a power source to charge the battery. A full battery will be indicated by a green LED. The supply voltage is 5V (USB-C). Estimated charging time is about 2 hours.

After the massage is completed, the device can be turned off by pressing the button again. The massager is ideal for people who require deep relaxation, want to improve the condition of their hair or improve the circulation of the scalp. With systematic use of the massage, you will reduce the symptoms of fatigue, insomnia or chronic headaches.

The strength of the massage should be adjusted to your own preferences. The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after August 13, 2005.

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The manufacturer hereby declares that the Medivon Neptun device complies with Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska  
Model: Medivon Neptun