medivon.

## User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK  
SL | EE | LT | LV | HU | RO | RU | NB | FI  
EL | HR | SV | DA | NL | TR | BG | AR

PL | Więcej informacji o produkcie dostępnych jest na stronie producenta pod adresem:  
EN | More information about the product is available on the manufacturer's website at:  
DE | Weitere Informationen zum Produkt finden Sie auf der Website des Herstellers unter:  
FR | Plus d'informations sur le produit sont disponibles sur le site Web du fabricant à l'adresse:

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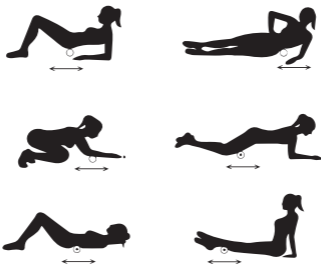
## Device operation


Before using it for the first time, make sure that the device is charged, then start the massager by pressing the POWER button, and then select the appropriate operating mode according to your preferences. When the massager is turned off, short press the button to turn it on. The default operating mode will start. Press the button to cyclically change the operating modes and massage intensity. The mode will be signaled by the LED.

The device has 3 operating modes (low - medium - high) and intermittent mode.

If the device is discharged, connect the device to a power source to recharge the battery. A full battery will be indicated by 3 lit LEDs. The supply voltage is 5V / 1A (USB-C).

The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use.

Below are examples of positions for using the massager to effectively massage individual muscle groups.

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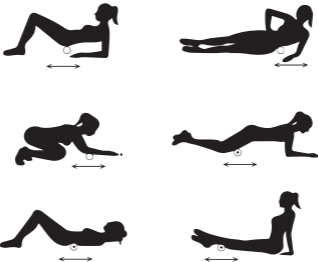Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. This device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the RecoRoll Trio device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A. | ul Grzybowska 87,  
Warszawa, 00-844, Polska  
Model: RecoRoll Trio