
medivon.

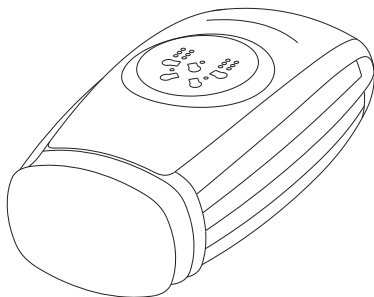Hand massager

Air massage technology

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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1. Operation of the device

The device is very easy to use - one press of a button is enough to start a relaxing hand massage which will greatly improve your blood circulation and relax tense hand muscles after a hard day's work.

Put your hand into the massager and then hold the POWER button for about 2 seconds to start. The device's default timer is 15 minutes - after this time, the device will turn off automatically. You can also turn them off by holding down the POWER button again.

Changing the massage intensity - to change the massage level, just press the appropriate button on the left.

The intensity will change cyclically from the weakest to the strongest L1 - L2 - L3 - L4 etc.

Changing the massage mode - the massager has defined 6 automatic massage modes. By pressing the MODE button, we can select cyclically changing programs, the one that currently suits us best.

The vibration option is intended for people who appreciate even better effects of relaxation and improved blood circulation in the hands. Using the vibration button, we can change their strength or turn it off completely.

Warm-up function - the built-in heads act as a warming compress, this function is activated by the heating button on the right side of the control panel. After starting the warm-up function - the temperature in the device will start to rise gradually.

More information about the product can be found on the manufacturer's website at:
www.medivon.pl/do-pobrania

2. Dangers

1. Do not use spare parts that are not approved by the manufacturer.
2. Unprofessional service should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid a hazard.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After completing the massage, disconnect the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, please stop using the massager immediately.
8. If the temperature is too high or the product starts to smell bad or smells burnt, stop using it immediately and disconnect the power.
9. Do not carry the product by pulling on the power cord to prevent the power cord from loosening.
10. This product is not intended for people who have no physical, sensory or mental abilities or have no experience and knowledge in this field (including children).
11. People insensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver.

3. Warnings

If you have any of the following conditions, consult your healthcare practitioner before use:

- Pregnant women, children or the elderly
- People who suffer from skin infections or have reddened and burned skin
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are treated or who experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

4. Cleaning and maintenance

Turn off the product and unplug the power cord before cleaning. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with harsh detergents such as benzene, thinners, and store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or product damage.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after 13.08.2005.

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The manufacturer hereby declares that the Medivon Hand device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, and Commission Delegated Directive (EU) 2015/863 of March 31, 2015 amending Annex II to Directive 2011/65/EU. The full content of the declaration of conformity is available at the following address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Hand